



Android

# MuseShape2

Multi-functional Body Shape &  
Skin Tightening Machine

Vacuum+RF+Roller+IR+Cavitation

# DISPLAY



# Treat yourself to a sexier figure!



Body Contouring



Cellulite Removal



Skin Tightening




Body slimming




# Human-Machine Interaction

Use Android smart operating system

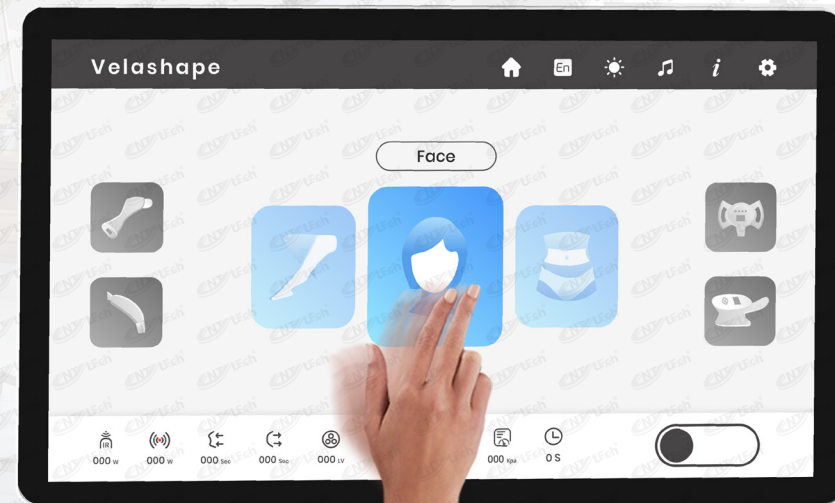

# ANDROID SOFTWARE



## Support system update

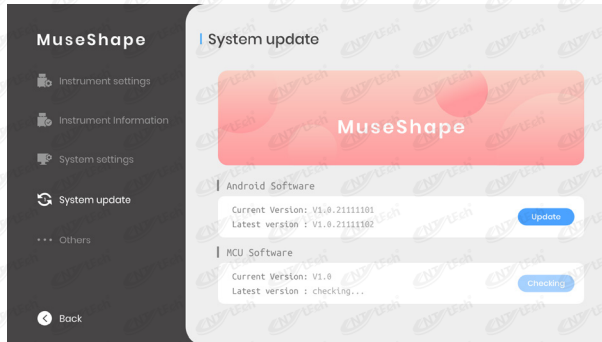

## Support Internet access



## Customer information storage and processing

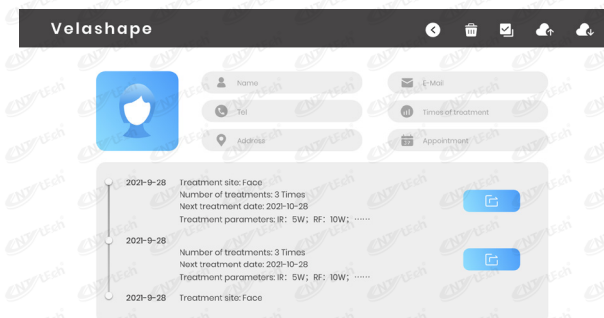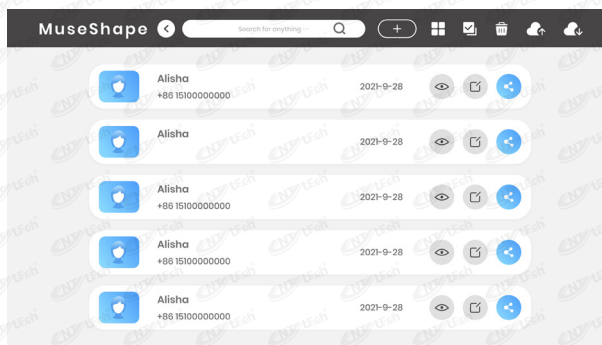

# HANDLE



## Handle 1

Treatment area:  
Abdomen.



## Handle 2

Treatment area:  
Back, Abdomen, Buttock.



## Handle 3

Treatment area:  
Chest, Back, Abdomen,  
Buttock, Leg, Arm.




## Handle 4

Treatment area:  
Face, Eye, Neck.




# HANDLE

Handle 1

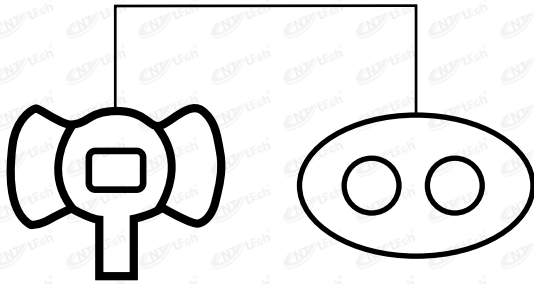


90mmX120mm

Handle 2

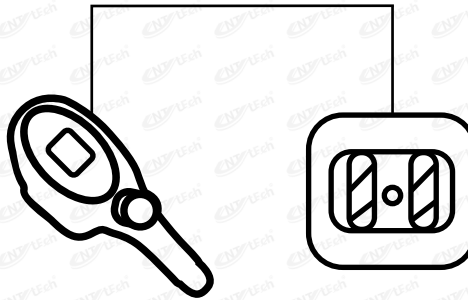

40mmX66mm

Handle 3

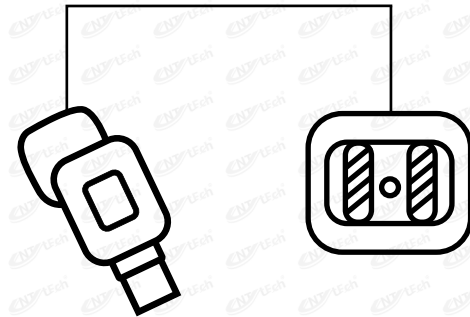

30mmX44mm

Handle 4



4.5mmX9mm

9mmX13mm

16mmX25mm

# HANDLE

## Touch screen

---

During treatment could adjust the parameters depending on customers' situation.



## Light weight

---

Big applicator weight is only 950g

Small applicator weight is only 350g

Reducing operator's hand pressure



# APPLICATION



# ADVANTAGE

five technologies in one machine: Vacuum+RF+Roller+IR+Cavitation.


4 unique and luxurious handles can treat whole body and face, over 20 projects supported.

## Advantage

Bi-Polar RF 1MHz/4MHz, different RF can be set for different body parts, more comfortable and better effect.

Perfect control system: the adsorption strength, adsorption frequency, roller speed, roller working direction, RF power. etc can be adjusted freely to meet the various treatment demands of patients.

# ADVANTAGE

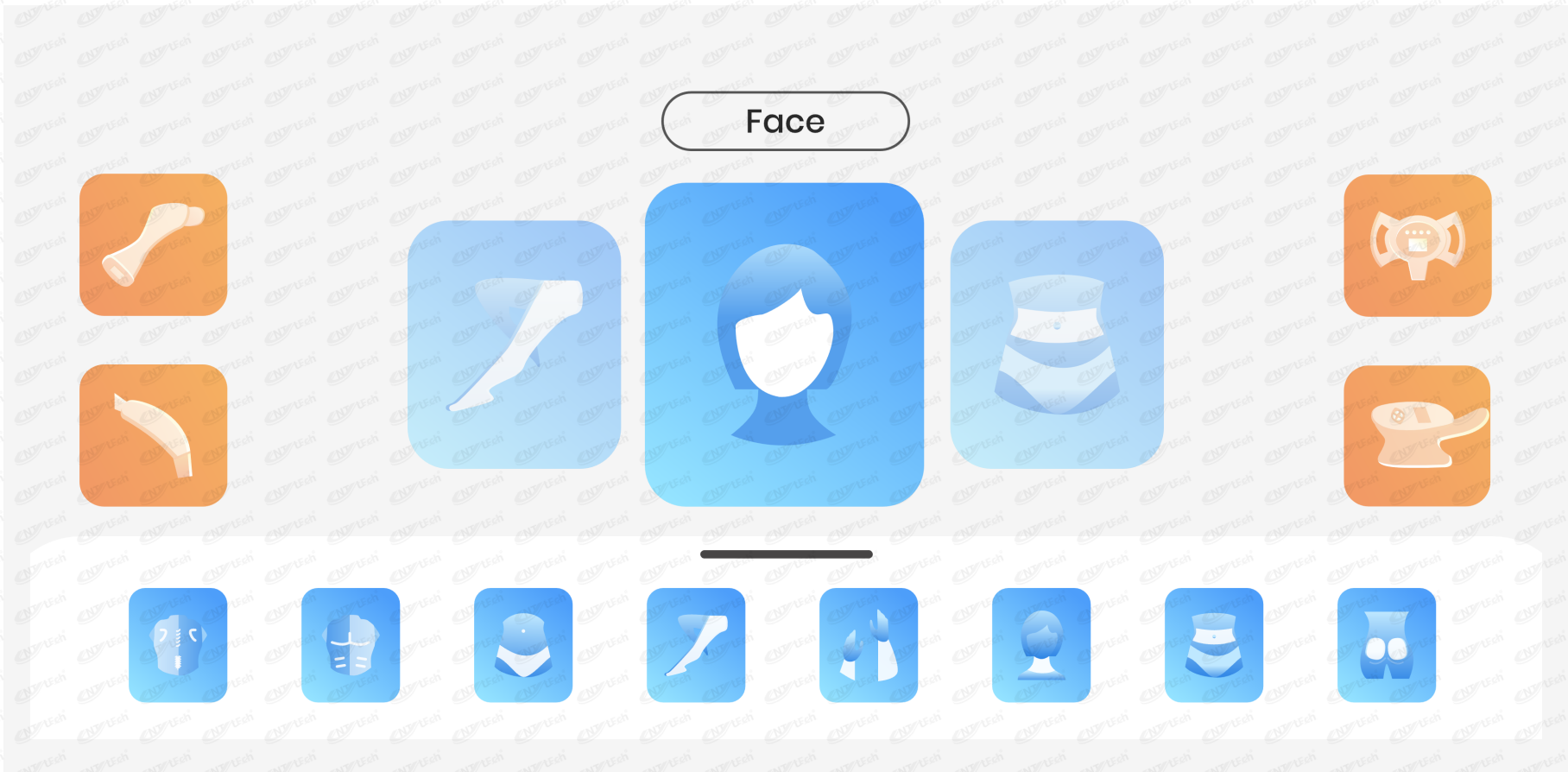

## Personalized Treatment

Select the appropriate applicator for different treatment parts  
Customized a personalized treatment plan.

**All body**

# ADVANTAGE

## Feel the details through heart

The Ultra silent system of MuseShape improve user comfort, easy to do the body shaping care.


The sound of flip book  
20-40dB

MuseShape running sound  
40-60dB

Speak loudly  
70-80dB




# ACCESSORIES

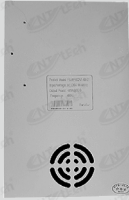

Power supply

Power supply: Support Bi-Polar RF 1MHz/4MHz, more comfortable RF treatment, better treatment effect.

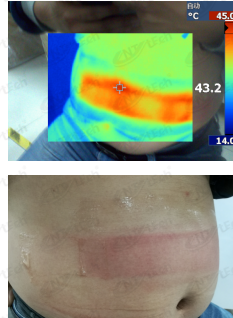

Japan imported ULVAC pump

Japan imported ULVAC pump: gas flow 66 kpa, no noise, longer life, high comfort of negative pressure adsorption.




Electromagnetic valve

South Korea imported YSC pneumatic Electromagnetic valve: Max 1.5Mpa, long life.



Proportional valve

Japan imported SMC proportional valve: real time pressure regulation.

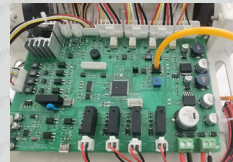

Imported pressure transducer


Imported pressure transducer: Small size, good temperature characteristics, high reliability.



940nm laser



940nm IR testing

# THEORY

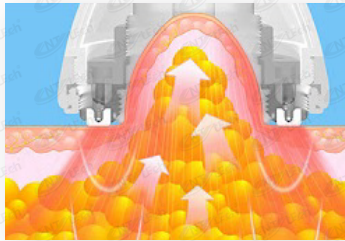



Vacuum—Regularly vacuum massage can promote lymph draining and decrease volume of fat cells.

RF—RF stimulates blood microcirculation, strengthen lymph system, promote fat Metabolism. Stimulation of thermal energy can reach to 0-15mm of skin, cause blood relax. increase body metabolism speed, enhance ligament and muscle flexible.

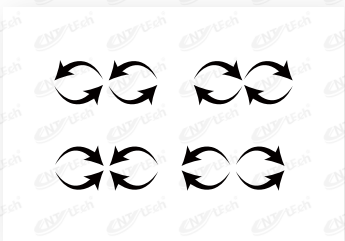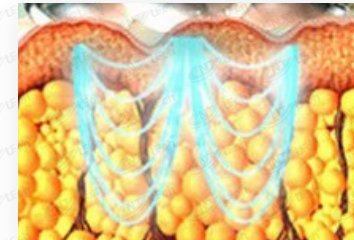

Roller - Roller massage of four modes, break up adipose connective tissue.

IR-940nm Infrared light Increase oxygen levels of body, make inert water transformed independent water molecules.

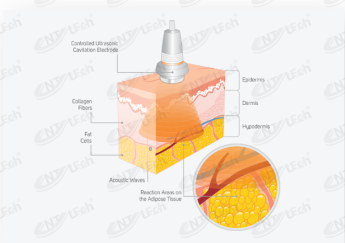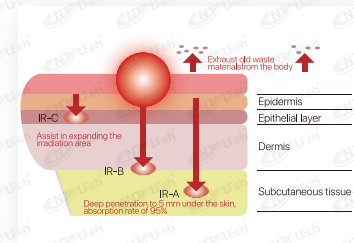

Cavitation-Ultrasonic cavitation, or ultrasound cavitation, is a cosmetic procedure that's used to break apart fat deposits in your body.

# THEORY

Infrared RF Vacuum Roller technology combines infrared light, bi-polar radio frequency energy and vacuum, which cause deep heating of the fat cells, their surrounding connective tissue and the underlying dermal collagen fibers. This type of efficient heating and vacuum stimulates the growth of new and better collagen and elastin which results in localized reduction in skin laxity, body volume, and an overall improvement in skin structure and texture.



Before



Affter



Before



Affter

# PARAMETER

Items	Specifications
LCD screen	1)Display screen: 13" Android intelligent screen 2) Display screen of handpiece No.1, No.2: 2.8" No.4: 3.5"
Safety checking	Real time on line
Negative pressure	1)Pressure regulating range: 1-8 level 2) Relative value: 30kPa -60kPa
Roller	1)Rotate speed: 1-5 level 2)Working mode: 4 types
RF	1)Frequency: 1MHz / 4 MHz 2)Energy density: Max: 60J/cm3
Ultrasonic	1)Frequency: 40kHz 2)Power: MAX60W
Number of handpiece	4
Treatment area	4.5mm*9mm, 9mm*13mm, 16mm*25mm; 30mmx44mm,40mmx66mm, 90mm*120mm
Mode of power supply	AC220V , 50Hz/AC110V , 60Hz
Machine size	420*400*1140mm
Package size	600*500*1180mm
N.W. /G.W.	60kg/75kg
Packing method	Aluminum case
Packing list	No. 1 big vacuum RF laser head for body * 1 No. 2 vacuum RF laser head for arms * 1 No. 3 vacuum RF laser head for face * 1 No. 4 dual vacuum cavitation head for body * 1 Power cord * 1





# COMPANY PROFILE



Ample spare parts reserve.  
Rich shell options.  
Instant delivery time.



UNT has a strong research and development team, including electronics engineers (designer; AVR programming; circuit layout) construction engineer, UI designer engineer QC control engineer...etc. We provide OEM/ODM design and program update.



In UNT we have been at the forefront of laser diode technology since 2010, We are committed with the clinics and practices to deliver them a platform that gives them the ability to get there results that the customers demand.



## Beijing UNT Technology CO.,Ltd

Address: M2-1 Area, Xinggu Development Zone, Pinggu District, Beijing 101200, China.

Website: <https://untlaser.en.alibaba.com/>

[www.untlaser.com](http://www.untlaser.com)

Contact: 86-010-61960002

